

1

Introduction – An Overview of Practical Electrochemistry

Students in natural sciences as well as professionals in numerous areas will meet electrochemical methods, concepts, and processes in many fields of science and technology. Accordingly, any conceivable selection of possible experiments intended as an illustration of this width, the numerous possibilities of electrochemistry, and an introduction to the subject has to be similarly broad. According to the book's purpose and intention this will be achieved by the width of the selection of experiments, the scope of the practical (instrumental) requirements, and the necessary level of knowledge. Convenient use of the book and logical arrangement of the essentials of the theoretical introduction suggest a rational arrangement of experiments. As already proposed and executed elsewhere (R. Holze: *Leitfaden der Elektrochemie*, Teubner, Stuttgart 1998 and *Elektrochemisches Praktikum*, Teubner, Stuttgart 2000), electrochemistry in equilibrium, i.e. without flow of current and conversion of matter, is followed by electrochemistry with flow of current. In the first chapter, measurements of electrode potential and their application in, e.g., the determination of thermodynamic data are treated. The second chapter deals with all kinds of experiments where an electrical current crosses the electrochemical interface. Applications of electrochemical methods (both without and with flow of current) are handled in a chapter on electrochemical methods of analytical chemistry¹⁾. This chapter also contains experiments helpful in elucidating, e.g., mechanisms of electrode processes (a somewhat broader meaning of analytical) if the focus of the experiment is not on the experimental method itself, thus suggesting its inclusion in one of the preceding chapters. According to the growing impact of non-traditional, in particular spectroscopic, methods in electrochemistry a small section of experiments from this branch follows; unfortunately, the feasibility of these experiments depends crucially on the presence of mostly expensive instrument. Electrochemical methods of energy conversion and storage are of utmost practical importance, and numerous first personal interactions with electrochemistry

1) The term analytical electrochemistry is neat and seemingly well-defined, but unfortunately in daily life its use is somewhat confusing: The term "analytical" is sometimes applied to qualify a certain branch of elec-

trochemistry (in contrast, e.g., with synthetic electrochemistry); sometimes it means application of electrochemical methods in analytical chemistry – as intended here.

deal with these devices. They combine applied aspects of equilibrium (i.e. thermodynamic) electrochemistry and electrochemical kinetics; consequently in this chapter only those experiments are collected where these aspects are not dominant. Electrochemical methods in industrial (synthetic) chemistry applied in the production of base chemicals like e.g. chlorine or sodium hydroxide and in synthetic procedures are subject of some experiments in the final chapter, the difficulties of the transfer of a large-scale industrial process into a simple laboratory experiment limit the selection.

The following descriptions of experiments are organized according to a general scheme. A brief statement of the experimenter's task and the aim of the experiment are followed by a condensed description of the theoretical foundations, essentially for understanding the experiment. This information cannot replace the respective parts of a textbook or the original reports in the primary literature. Besides references to primary sources, the respective sections of C. H. Hamann, A. Hamnett, and W. Vielstich: *Electrochemistry* (Wiley-VCH, Weinheim 2007) are quoted as EC and the respective page number: EC:xx. Some methods like polarography and cyclic voltammetry are employed in several experiments; nevertheless their fundamentals are described only once when the method is introduced first. No attempt is made in the descriptions to list all conceivable applications of the method used in this experiment. The tempting concept to arrange experiments according to difficulty or complexity of the experimental apparatus was discarded soon as being too personal and subjective. Instead, the reader and user of this book will easily select experiments according to his or her personal interests and intentions; the comparison of available and necessary equipment can subsequently be performed easily as well as the estimate of the required knowledge for successful execution.

The description of the execution of an experiment starts with a list of necessary instruments and chemicals. Possible alternative instrumentations are highlighted; the subsequent description is nevertheless limited to one experimental way only. The description contains – if necessary – a schematic circuit diagram of the setup and sketches of the construction of the apparatus or parts of it. The execution of the suggested measurements is briefly outlined. Potential pitfalls and unusual details are indicated. The way from the raw data to the desired results is sketched. Final questions including those pertaining to the practical execution help to confirm the newly acquired knowledge. Extensive calculation training examples are not included, these can be found in the textbook by J. O'M. Bockris and R. A. Fredlein. Typical results are displayed without cosmetic tidying up, this will encourage the user, it also demonstrates the level of skill needed to obtain satisfactory agreement between literature data (always quoted according to bibliographic standards from generally available textbooks for comparison) and one's own results.

Practical Hints

In most experiments aqueous solutions are used. If not stated otherwise ultra-pure water (sometimes called 18 M Ω -water because of the typical specific resistance value of this water) is used. It can be obtained by afterpurification of deionized water by various commercially available purification systems. As an alternative doubly distilled (bidistilled) water can be used. In some experiments simple deionized water can be used. Because especially in demanding conductivity²⁾ and potentiometric measurements traces of impurities present in deionized water may cause erroneous results, blind tests are required, in particular when water of less than ideal purity is used. In some cases not only the desired concentration of a necessary solution but also the amount of the selected chemicals needed for preparing the requested amount of solution are given for ease of preparation. When cells or other experimental setups with volumes different from the suggested setup are used these numbers must obviously be corrected. Purification of organic solvents has been thoroughly described by Mann (C.K. Mann: *Nonaqueous Solvents for Electrochemical Use*, *Electroanalytical Chemistry 3* (A.J. Bard, ed.), Marcel Dekker, New York 1969, p. 57); further information on electrolyte solutions based on organic solvents has been collected by Gores and Barthel (H.J. Gores and J.M.G. Barthel, *Pure & Appl. Chem.* **67** (1995) 919).

Electrodes

As suggested by W. Nernst, the term electrode should always be applied to a specific combination of an electronically conducting material (e.g. a metal, graphite, a semiconductor) and an ionically conducting phase in contact with this material (e.g., an aqueous solution, a polymeric electrolyte, a molten salt). The need for this use will be neatly illustrated in experiment 3.13 with lead being in contact with various electrolyte solutions; quite obviously the term lead electrode becomes ambiguous. In daily life the term electrode mostly refers to the electronically conducting component only. This well-established usage will not be completely suppressed in this book; nevertheless the possible confusion will be addressed repeatedly.

2) Meaning and use of the terms conductance and conductivity and their respective reciprocal relatives resistance and resistivity are slightly confusing and frequently inconsistent. According to standard reference works conductance applies to the conducting power of matter without requiring specific dimensions of the sample whereas conductivity refers to the conductance of a sample

with unit dimensions, i.e. it refers to specific conductance. The frequently employed term "specific conductivity" might thus be popular, but it remains a pleonasm or tautology. In this book the standard terms conductivity meter or conductivity measurements are employed – although in most cases only conductances are measured.

In some experiments electrodes of special shape and construction prepared from selected materials are needed; details are provided in the descriptions of the experiments. In many experiments electrodes of a fairly general type and construction will be used. Because they can be prepared easily in a glassblower's shop or even without any expert help, some suggestions are given below.

Frequently, metal sheets (of noble metals like platinum or gold) are used as working and counter electrodes. These electrodes can be manufactured easily by spot welding a metal wire to a piece of sheet metal (about 0.1 to 0.2 mm thick). After extending the metal wire with a piece of copper wire connected with hard solder (i.e. silver solder; soft solder is not recommended because it will most likely melt in the subsequent glassblowing operation; in addition soft solder forms alloys with gold making a reliable connection impossible) the noble metal wire is sealed into a glass tube. Glass with low melting point is preferred because the low viscosity of the molten glass obtained even at moderate temperatures provides a tight glass-metal seal. With platinum, borosilicate glass of higher melting point can be used. Particularly useful even for silver wires are lead dioxide-based glasses. The glassblower must avoid reducing conditions when operating his blowtorch. Unfortunately these glasses are hard to get. When a spot welding machine is not available simple metal wire spirals can be used instead of sheet electrodes. Instead, with a glass-metal seal the wire can also be fixed with epoxy glue; unfortunately this connection is mechanically and chemically less stable, may not be exposed to some aggressive cleaning solutions, and in addition may release traces of contaminants into the electrolyte solution.

In the case of a very simple preparation of a microelectrode this epoxy glue is essential as shown below (Fig. 1.1). A glass tube is heated until its end almost collapses and a very narrow opening remains. A carbon fiber is fed through the opening. The opening is closed with epoxy. By pulling and pushing the fiber gently it is coated with epoxy and completely surrounded with the glue. After curing the epoxy an electric connection to a copper wire is provided on the inside with conductive silver or graphite cement.

As reference electrode mostly hydrogen or metal ion electrodes are used. Saturated calomel and silver chloride electrodes (EC:99) are particularly popular. Figure 1.2 shows typical constructions.

Various aqueous as well as nonaqueous electrolyte solutions can be used; in all cases the high chloride concentration may result in contamination of the electrolyte solution in the electrochemical cell because of halide diffusion

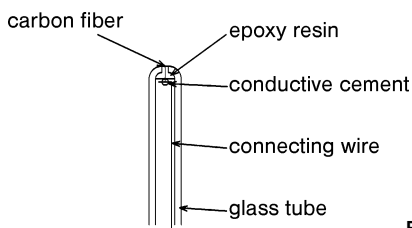


Fig. 1.1 Cross section of a simple microelectrode.

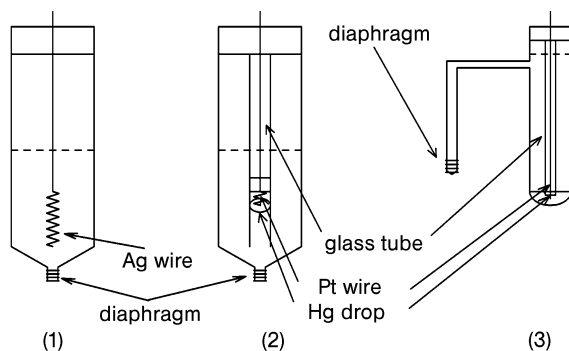


Fig. 1.2 Cross sections of various types of reference electrodes: (1) Silver/silver chloride electrodes; (2, 3) Calomel or mercurous sulphate electrode, dotted lines indicate solution filling level.

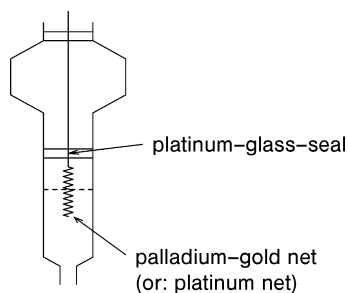


Fig. 1.3 Hydrogen electrode according to Will.

through the separating diaphragms etc. In alkaline solution disproportionation of calomel may occur.

A particularly simple design of a hydrogen electrode³⁾ as suggested by Will (F.G. Will and H.J. Hess, *J. Electrochem. Soc.* **120** (1973) 1; F.G. Will and H.J. Hess, *J. Electrochem. Soc.* **133** (1986) 454) is shown below (Fig. 1.3).

Platinum wire net (better: a palladium or palladium/gold wire net) is spot-welded to a platinum wire. This assembly is fixed inside the glass tube drawn into a capillary with a glass-metal seal. The active surface area of the net is enlarged by electrolytic coating with platinum black. Using the electrolyte solution of the planned experiment, an auxiliary electrode (e.g. platinum wire) as anode

3) Sometimes hydrogen electrodes are called reversible or relative hydrogen electrodes RHE. The former term is obviously redundant. The term reversible implies the presence of a reversible reaction: An electrode reaction proceeds on the same reaction pathway in both directions at significant rates; without this no stable electrode potential would be established. Alternatively the

term might imply reversible in a thermodynamic sense – being at equilibrium. Because a reference electrode is always used without any flow of current it is in equilibrium. The term relative refers to the fact that a nonstandard proton concentration may be present, thus the electrode is not a standard electrode.

and a DC power supply (a few volts will be enough), the metal net inside the capillary is charged with hydrogen by cathodic reduction of protons until a hydrogen gas bubble is formed. The bubble may stay in place for several weeks when the glass-metal seal is sufficiently tight; accordingly the electrode potential remains stable. In the case of organic compounds dissolved in the electrolyte solution, which may react at the platinum surface, potential shifts may occur as a consequence of electrode poisoning; in the case of neutral electrolyte solution the potential is also less stable because the exchange current density of the hydrogen electrode reaction is pH-dependent (EC:344). The activity of protons in the solution filled into the hydrogen electrode will be unity only in rare cases. Consequently this electrode may not be called a standard hydrogen electrode SHE (the term normal hydrogen electrode NHE should be avoided altogether because the term normal may be taken as designation of a certain concentration); because the hydrogen electrode potential is related to the proton activity this hydrogen electrode is sometimes called relative hydrogen electrode RHE.

With nonaqueous electrolyte solutions the use of reference electrodes containing aqueous electrolyte solutions is still possible with salt bridges (see below) providing the ionically conducting connection. The diffusion potentials created at the phase boundaries between the various solutions may cause considerable experimental errors if not corrected properly. In addition the slow diffusion of solution components (including of course water) into the nonaqueous electrolyte solution may cause undesirable chemical reactions or other experimental artifacts. Thus electrolyte solutions employing nonaqueous solvents are an attractive option. Unfortunately electrodes of the second kind (e.g., silver/silver chloride, for details see above) are prone to disproportionation and finally decomposition resulting in reference electrode potential drift. In addition the potential of these electrodes depends on the employed solvent and thus comparison between experimental results obtained with different solvents may be unreliable. Redox systems with a formal potential E_0 defined as the midpoint potential between their respective oxidation and reduction peak potentials as observed in a cyclic voltammogram ($E_0 = E_{p,red} + (E_{p,ox} - E_{p,red})/2$, for details see experiment 3.18) have been employed as point of reference repeatedly provided that the value of E_0 does not depend on the solvent. Ferrocene (see experiment 3.18) has been suggested as a candidate (R. R. Gagné, C. A. Koval, and G. C. Lisensky, *Inorg. Chem.* **19** (1980) 2854) because the iron ion in the center between the two cyclopentadienyl ligands seems to be shielded fairly well from the solvent. Unfortunately it has turned out that upon oxidation of ferrocene most of the charge is removed from the ligand; the actual charge on the iron is changed only by about $1/10^{\text{th}}$ of an electron charge. As a substitute decamethyl ferrocene or decaphenyl ferrocene have been suggested; for a review see I. Noviadri, K. N. Brown, D. S. Fleming, P. T. Gulyas, P. A. Lay, A. F. Masters, and L. Phillips, *J. Phys. Chem. B* **103** (1999) 6713. In an actual experiment a reference electrode of the second kind filled with a nonaqueous electrolyte solution may be used; at the end of the experiment some ferrocene is added, a cyclic voltammogram is recorded, and all electrode potentials are converted to this reference (e.g. the fer-

rocene) scale. Only in this way can results obtained with different solvents at different places be made compatible and comparable and comparisons be made possible.

Although a large selection of electronic reference voltage sources with a precision sufficient for calibration purposes is available (for a suggested circuit see the appendix) electrochemical reference cells are still in use. The only infrequently encountered Clark cell ($\text{Zn}|\text{ZnSO}_{4,(\text{aq.}, \text{sat.})}|\text{ZnSO}_{4, \text{sol.}}+\text{Hg}_2\text{SO}_{4, \text{sol.}}|\text{Hg}$) has been replaced by the Weston cell ($\text{Cd}|\text{CdSO}_{4,(\text{aq.}, \text{sat.})}|\text{CdSO}_{4, \text{sol.}}+\text{Hg}_2\text{SO}_{4, \text{sol.}}|\text{Hg}$). The cell reaction of the latter has significantly lower entropy of reaction resulting in a smaller temperature coefficient of the cell voltage.

Measuring Instruments⁴⁾

Aside from special instruments typical of certain experimental methods (e.g., polarographs for polarography), simple electronic instruments are frequently needed in particular for the measurement of voltage and current. Standard analog or digital multimeters are sufficient in most cases. During measurements of electrode potentials performed as the measurement of the voltage between a reference electrode and the working electrode under study the flow of electric current must be avoided. Ideally this has been done using a compensation circuit. This procedure is cumbersome and today of only limited practical importance. Voltmeters with very high input resistance ($R_i > 10^{12} \Omega$) are sufficiently close to this ideal method. During selection of instruments for precision measurements this input resistance merits particular attention. Budget-priced multimeters are frequently built with a voltage divider circuit for range selection at the input resulting in fairly low input resistance values; these should be used with appropriate respect. With commercially available digital voltmeter modules (no input voltage divider added) with a 2-V range a powerful instrument for potential measurements can be built at low cost. Measurements of current-potential curves require current meters with a very small measuring (shunt) resistor; ideally its value should be zero. This can be realized with simple circuits based on operational amplifiers (current followers) which are particularly helpful during measurements with electrochemical cells providing only small output voltages (fuel cells, batteries). Because of lower price and apparently higher precision, digital instruments are frequently preferred. During experiments wherein a voltage has to be adjusted to a certain value an analog instrument may be better because trends in the measured signal can be discerned more clearly. Digital instruments with an additional bargraph display may be a compromise; it may be difficult to become accustomed to the flickering bargraph display.

4) A complete setup including a potentiostat built as plugin for a standard desktop PC, a cell, several electrodes, and software suitable for running several of the experiments

described below including a workbook is offered by Sycopel Scientific Instruments, 15 Sedling Road, Washington NE38 9BZ, Great Britain.

Electrochemical Cells

As well as specific cells designed just for a certain experiment as described in the text related to it, some generally used types of electrochemical cells have come into use. A simple beaker will be sufficient only in a few cases because operation with inert gas atmosphere above the electrolyte solution and sufficient purging of the solution itself are impossible with this arrangement. In addition reliable and reproducible mounting of electrodes can be done only with additional holders etc. For cyclic voltammetry a cell depicted below in cross section, commonly called H-cell (because of its shape) is frequently used (Fig. 1.4).

For electrochemical impedance measurements and further experiments with AC-modulation of the electrode potential spherical working electrodes or circular disc-shaped electrodes embedded in an inert material mounted in the center of a symmetric cell vessel are advantageous. Additional ground glass feed-throughs hold glass tubes for gas purge of the electrolyte solution and the counter and reference electrode. Reduced exchange of electrolyte solution between the interior of the latter tubes and the main compartment of the cell can be achieved by closing the tubes with porous glass or ceramic frits. This type of cell can also be employed in studies with a rotating disc electrode as shown below (Fig. 1.5).

In precision studies⁵⁾ and in experiments, where mixing of electrolyte solutions must be prevented, “salt bridges” (Fig. 1.6) are needed. They provide electrolytic connection between cell compartments without bringing the electrolyte solutions in these compartments into direct contact. In a very simple design a piece of plastic tube filled with a suitable electrolyte solution and closed with cotton wool plugs or plugs made of filter paper may be sufficient. The electrolyte could be 1 M KNO_3 , both ions having similar mobility (P.W. Atkins, J. de Paula, *Physical Chemistry*, 8th ed., Oxford University Press, Oxford 2006, p. 1019) which results in an approximate compensation of diffusion potentials (P.W. Atkins, J. de Paula, *Physical Chemistry*, 8th ed., Oxford University Press, Oxford 2006, p. 216; EC:112;146) generated at the interfaces the cell solutions.

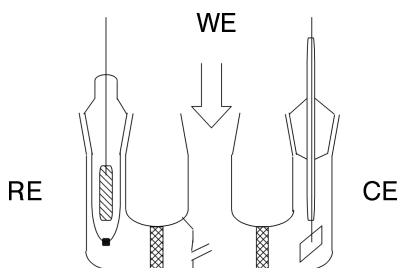


Fig. 1.4 Cross section of an H-cell for electrochemical experiments, WE: working electrode, RE: reference electrode, CE: counter electrode.

5) For example measurements of standard electrode potentials or cell voltages wherein diffusion potentials must be either avoided or precisely known.

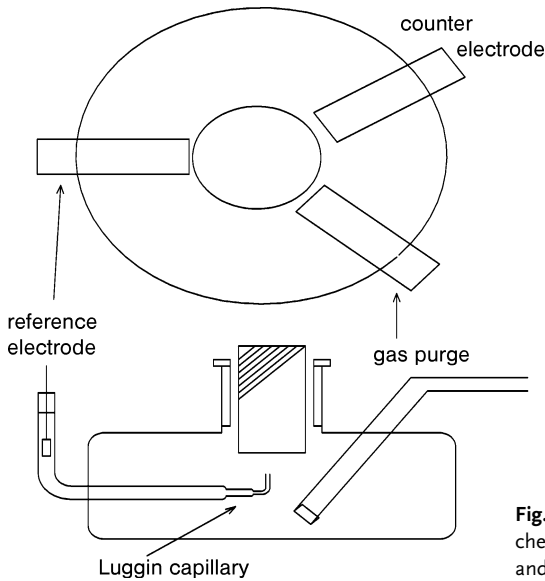


Fig. 1.5 Cross section of an electrochemical cell for AC-measurements and rotating disc electrode studies.

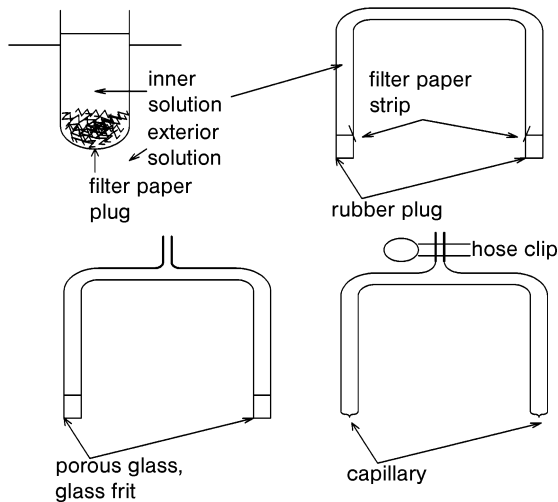


Fig. 1.6 Cross sections of some salt bridge constructions.

In addition both ions show only minor interference with most electrochemical processes and cell reactions. More reliable and stable are ground glass diaphragms or porous glass plugs (made of "VYCOR[®]"⁶). Some typical designs are shown above (Fig. 1.6).

⁶ These plugs must be kept wet, as drying might result in destruction.

Data Recording

Extensive data recording may not be necessary with most experiments. Only in studies of processes where data have to be recorded at a high rate (e.g., in cyclic voltammetry or spectroscopy) are recording devices needed. In addition to oscilloscopes (which provide printed output in most cases only with added hard- and software at a considerable price), analog recorders (X-Y-recorders or X-t-recorders) have been used. The omnipresence of cheap and powerful computers and the steep rise in the price of analog recorders suggest the use of computerized data recording setups. Easily programmable and cheap ADDA-interface⁷⁾ cards enable the generation of input signals for electrochemical experiments (via the DA-capability) and the recording (via the AD-capability) making the computer a replacement for both recorder and function generator. Subsequent data storage and treatment would preferably be by the use of a computerized workplace instead of cumbersome paper shuffling. Some drawbacks should be kept in mind nevertheless.

Commercial software for electrochemical experiments is powerful and expensive in most cases (for an exception see footnote on p. 8), so that its purchase just for a laboratory course may be inappropriate. In addition the use of analog function generator and X-Y-recorder has some advantage when adjusting experimental parameters like, e.g., potential limits, scan rates etc. because the effect will become visible immediately. With computerized setups this is impossible in most cases; thus the direct cause-effect relationship important in education and hands-on teaching gets lost. Interface cards need careful hard- or software calibration, which tends to be forgotten frequently. Because calibration needs external voltage sources providing unusual voltages (not the values of standard cells described above) a simple circuit based on a high-precision voltage reference circuit is provided in the appendix. This circuit needs only a single initial calibration versus a high-precision reference. When selecting measuring ranges on potentiostats and other devices the input voltage range of the DA-interface card must be kept in mind. The range is frequently 2 V or 5 V. In potential measurements (actually the voltage between the working and the reference electrode is measured by the potentiostat) this range is a very good match. In the case of current measurements, frequently only small currents are observed. Unsuitably small shunt resistors provide only very small voltages resulting in a poor signal-to-noise ratio on the current display. Thus a shunt resistor as high as possible should be selected; in modern potentiostats the built-in current follower permits this good match.

7) ADDA: analog-digital/digital-analog converter